Product information

Extrusion set

2 magnets and 2 position discs

Maximum force: ca. 1 N


Magnet diameter: 3.80 mm

Magnet height: 2.45 mm

REF: S.62.01.Y1

Extrusion magnet for the **template**

REF: V.62.01.Y245.C

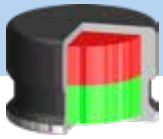Position discs (Spacer)


Height: 1.0 mm (0.33 N)
REF: P.62.01.Y100

Height: 2.0 mm (0.13 N)
REF: P.62.01.Y200

Extrusion magnet for the **root**

REF: V.62.01.Y245.R

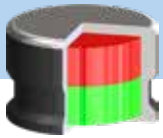
Further products

Titanmagnetics®

Magnetic attachments for implants:

- secure prosthesis hold
- easy handling
- no matrix change

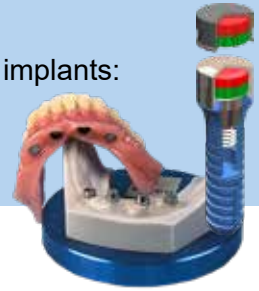Root cap magnets

For moulded root caps:
Root residues simply and quickly as an additional retaining element

StecoGuide

Precise implant planning:
Drill sleeves for every case

Please contact us at info@steco.de

steco-system-technik GmbH & Co. KG
Kollaustr. 6 • 22529 Hamburg • Germany
Tel. +49 40 55 77 81-0• www.steco.de/en

K.62.00.EN02 (03.2025)

steco system-technik

Extrusion magnets

The magnetic extrusion.
Simply attractive!

Titanmagnetics®

Problem

Tooth broken off?

Magnetic force lifts root and gingiva

Principle

The magnetic principle

- one magnet is fixed on the root
- another magnet is fixed on a template or a provisionary denture in a small distance
- the magnetic attraction between both magnets pulls the root out of the alveole slowly until both magnets meet
- bone tissue below the root is build up
- the soft tissue follows the root coronally
- the maximum force between the magnets is about 1 N (100 g)

Specific Equipment


The distance between the magnets is defined by magnetic position discs.

Solution

Crown instead of extraction

The root is factured too low to be used for crown fixation?


With magnetic extrusion the root can be pulled out of the alveole to be used as a crown base.

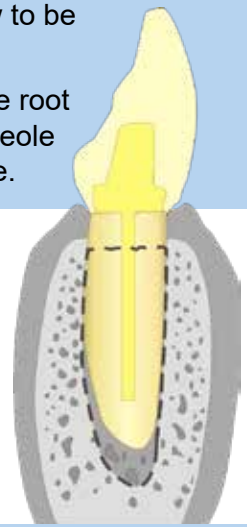Increase bone before implantation

With the help of magnetic extrusion patients bone and gingival tissue can be pulled upwards.

The implant can be inserted into more body‘s own bone tissue and the prosthetic solution profits from higher gingival level.

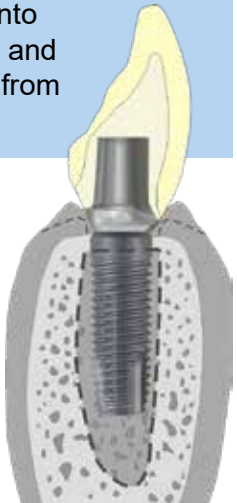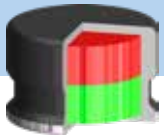
For more information, visit www.steco.de/en

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